


# rayence



10x12 Wireless Flat Panel Detector

## 1012WCC

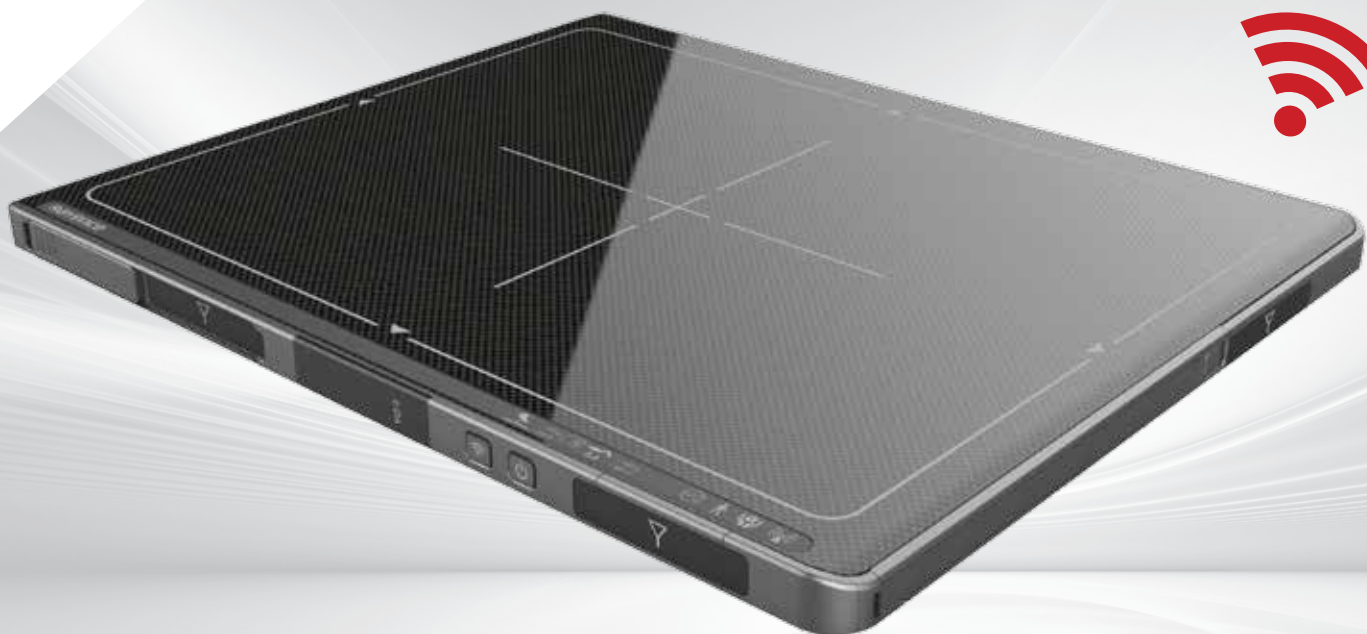

The ergonomically designed Rayence C-Series Cesium Iodide wireless detectors are designed to offer high levels of handling, functionality, and exceptional diagnostic image quality in the X-ray room and beyond

C  
series

## FEATURES

- ▲ Tapered, Recessed Edges
- ▲ Thin and Lightweight
- ▲ IPX6 Water Resistant
- ▲ Image Storage: 200 Images
- ▲ Room Sharing Functionality
- ▲ Auto Triggering Technology



### Superb Image Quality

1012WCC high Detector Quantum Efficiency (DQE) achieves superb image quality with low patient dose.



### Lightweight & Fast

1012WCC weighs only 3.7lb. Image preview occurs in less than 2 seconds.



### High Visibility OLED

Illuminated OLED window brightly indicates flat panel detector status to the user.



### Ergonomic Design

Curved edges and a non-slip surface makes lifting and handling easier for all tabletop extremity exams as well as for pediatric and neonate imaging.



### Durability

Supporting up to 660 lb., the 1012WCC is manufactured with a seamless magnesium, unibody construction and is combined with a shock, vibration, and scratch resistant carbon fiber composition.



### Water Resistant (IPX6)

1012WCC is water resistant to most typical water spills in a hospital as well as in outdoor applications.


**“ Providing patient throughput in your hospital and beyond ”**




### 1012WCC/WGC Specifications

ScintillatorType	1012WCC : CsI:Tl	-
Dimension	12.9 x 10.6 x 0.6	in
Weight	3.70 (incl. battery)	lb
Total Pixel Area	11.5 x 8.9	in
Pixel Pitch	127	µm
Effective Pixel Matrix	2264 X 1752	Pixels
A/D Conversion	14 / 16	bits
Preview time	≤2 (2x2 binning)	sec
Energy range	40 ~ 150	kVp
Pressure	Distributed : 661 Point : 330	lb
Limiting Resolution	Min. 2.5 / Max. 3.93	lp / mm
Battery OperatingTime	Typ. 4	Hours

### Detector Room Sharing

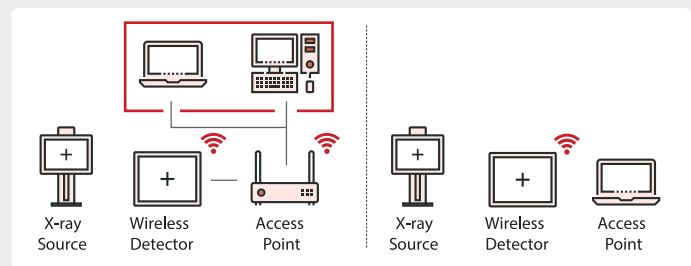

### Options

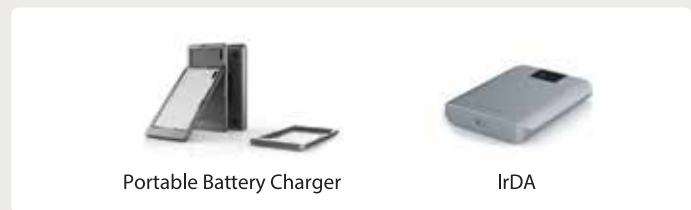

## Contact us

Rayence Inc.  
 81 Ruckman Rd, Unit A,B  
 Closter, NJ 07624  
 Office: 201.585.0290 Fax: 201.585.0293  
 Email: information@rayenceusa.com  
 www.rayenceusa.com

**rayence**  
 1012WCC